

AMREX[®]

AFFORDABLE ELECTROTHERAPY EQUIPMENT

Z-STIM 100[®] T.E.N.S.

Quality, Effective Therapy for Pain Relief



AMREX Z-STIM 100 T.E.N.S. Therapy

Transcutaneous Electrical Nerve Stimulation (T.E.N.S.) modality provides a sensory stimulation intended for symptomatic relief of a large number of painful syndromes. Prior to the cause of the pain being isolated, the T.E.N.S. modality relieves chronic intractable pain syndromes and eases pain where analgesic drugs would be contraindicated.

The Z-Stim 100 effectively relieves pain by applying stimulation through the self adhesive pad electrodes. The Z-Stim 100 can be placed in a shirt pocket, clipped to a belt or to a waistband. The Pulse Rate and Intensity are adjustable for optimal relief.

AMREX T.E.N.S. Application... The Affordable Alternative

The Z-Stim 100 fits your patient's budget as comfortably as it does their pocket. The Z-Stim 100 is backed by a two year parts and labor warranty and a tradition of over 60 years of American made electrotherapy experience.



Electrodes furnished may vary from those pictured

*"As simple to use
and convenient
to wear
as a pager"*



- Effective Therapy
- Quality Treatment
- Reliable Operation
- Durable Construction
- Easy to Use
- Affordable
- Wearable
- Concealable

800-221-9069

Z-STIM 100® T.E.N.S.


SPECIFICATIONS

Input Power Requirements

- Power Source 9 Vdc (standard 9 V battery)
- Current 20 mA

Electrical Stimulator

- Waveform Asymmetrical alternating current (no load)

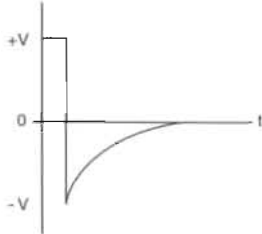

- Output Voltage 110 V peak (no load)
- Output Pulse Width 300 μ s at 50% V maximum (no load)
- Output Frequency 1 to 58 Hz

ACCESSORIES INCLUDED: 9 volt battery, patient leads, single patient reusable electrodes, belt clip, user's guide and storage case.

WARRANTY: AMREX provides a two (2) year warranty on the instrument to be free from defects in materials and workmanship under normal use and service.

FEATURES

- Two pad application
- Adjustable pulse rate
- Adjustable intensity
- Reusable pad electrodes included
- 9 volt battery included

CONTRAINDICATIONS AND WARNINGS

1. Do not stimulate over the carotid sinus nerves.
2. Do not use on patients with demand type cardiac pacemakers.
3. Do not apply transcerebrally.
4. Do not apply until the pain syndrome is diagnosed and the etiology is established.
5. The safety of TENS devices for use during pregnancy has not been established.
6. TENS therapy should only be used under the continued supervision of a licensed health professional.
7. TENS is a symptomatic treatment and as such suppresses the sensation of pain which would otherwise serve as a protective mechanism on the outcome of the clinical process.
8. Keep out of the reach of children.
9. Electronic monitoring equipment, such as EKG monitors or alarms, may not function properly when stimulation is in use.

PRECAUTIONS

Sensitive individuals may receive isolated cases of skin irritation or burning at the site of electrode placement following long term application.

Caution: Federal (U.S.A.) law restricts this device to sale by or on the order of a licensed practitioner.

Other AMREX Quality Products

- Interferential Therapy Equipment
- Muscle Stimulators
- Ultrasound Units
- Ultrasound Gel
- Conductive Spray
- Flextrode® Electrode System
- Mobile Utility Cart
- Soft Carrying Case

Distributed by:

Nobody Knows Amrex Better Than Amrex



For customer product support, contact Amrex Technical Services, the factory product support facility.

ATS
Amrex Technical Services

Amrex Technical Services is the factory site for product information and service. Keep your Amrex equipment operating factory new and in compliance with FDA regulations. Amrex Service Technicians will perform electrical safety tests, mechanical inspections and upgrade your equipment and accessories. . . typically within 72 hours!

For prompt, competitively priced factory service give Amrex Technical Services a call. A knowledgeable customer service specialist can be reached at 800-221-9069.

7034 Jackson Street, Paramount, California 90723 USA
Phone: 310-527-6868 Fax: 310-365-7343
Web Site: www.amrexusa.com
E-mail: amrex@amrex-zetron.com

AMREX®
AFFORDABLE ELECTROTHERAPY EQUIPMENT

800-221-9069

Revised 8/12